

Unistik®3 is available from your GP on prescription or buy direct at www.medicalshop.co.uk or call 0800 731 6959.

Unistik® 3 Junior

COMING SOON ON PRESCRIPTION



NEW

Gauge: 28G Depth: 1.4mm
Average Blood Volume: 6µL
Low Blood flow

Box of 100: AT1022
Box of 200: AT1024

PIP: 367-0098
PIP: 367-0015

Unistik® 3 Comfort

AVAILABLE ON PRESCRIPTION



Gauge: 28G Depth: 1.8mm
Blood Volume Range: 5µL – 30µL
Low Blood flow

Box of 50: AT1045
Box of 100: AT1042
Box of 200: AT1044

PIP: 347-4020
PIP: 317-2277
PIP: 317-2285

Unistik® 3 Normal

AVAILABLE ON PRESCRIPTION



Gauge: 23G Depth: 1.8mm
Blood Volume Range: 30µL – 75µL
Medium Blood flow

Box of 50: AT1008
Box of 100: AT1002
Box of 200: AT1004

PIP: 347-4038
PIP: 317-2301
PIP: 317-2293

Unistik® 3 Extra

AVAILABLE ON PRESCRIPTION



Gauge: 21G Depth: 2.0mm
Blood Volume Range: 75µL – 125µL
High Blood flow

Box of 100: AT1012
Box of 200: AT1014

PIP: 317-2269
PIP: 317-2251

Capillary Blood Sampling

Unistik® 3



medical shop

making life easier



To buy direct or for more information call us on 0800 731 6959 or visit us at www.medicalshop.co.uk



OWEN MUMFORD
Making a World of Difference

UK / IRELAND SALES & MARKETING

Owen Mumford Ltd
Primsdown Industrial Estate
Worcester Road, Chipping Norton
Oxfordshire OX7 5XP, United Kingdom
T : 01608 645555 F : 01608 646546
E : info@owenmumford.co.uk

www.owenmumford.com

OUR ACCREDITATIONS INCLUDE:

- ISO13485 Certification
- Medical Device Directive 93/42/EEC
- CMDCAS 13485 (Canadian Health Authority)
- Environmental Management Systems ISO14001



AT/1042-LT/MCC01128



Making blood sampling less of a pain

Unistik® 3

Single use Safety Lancets

Featuring our unique Comfort Zone Technology®

AVAILABLE ON PRESCRIPTION



OWEN MUMFORD

Unistik® 3

Single use Safety Lancets

Simple, comfortable blood sampling available on prescription.

Unistik® 3 offers a convenient way to test your blood glucose levels. A simple one-click action activates the lancet and our unique **Comfort Zone Technology®** helps mask the pain.

Safe and Hygienic.

Unistik® 3 is a single use lancet that cannot be reused reducing your risk of infection as you use a new device for each blood sample. The needle is also always shielded so there is no chance of an accidental injury either before or after sampling.

A lancet for everyone.

Unistik® 3 comes in four different variants enabling you to use the device which is right for your skin type.




Key benefits:

- SAFE – needle retracts immediately after use
- STERILE – single use, sterile device
- ACCURATE – targeted positioning
- LOCKOUT INDICATOR – prevents re-use

Sterile.

The twist-off protective cap keeps the lancet sterile.

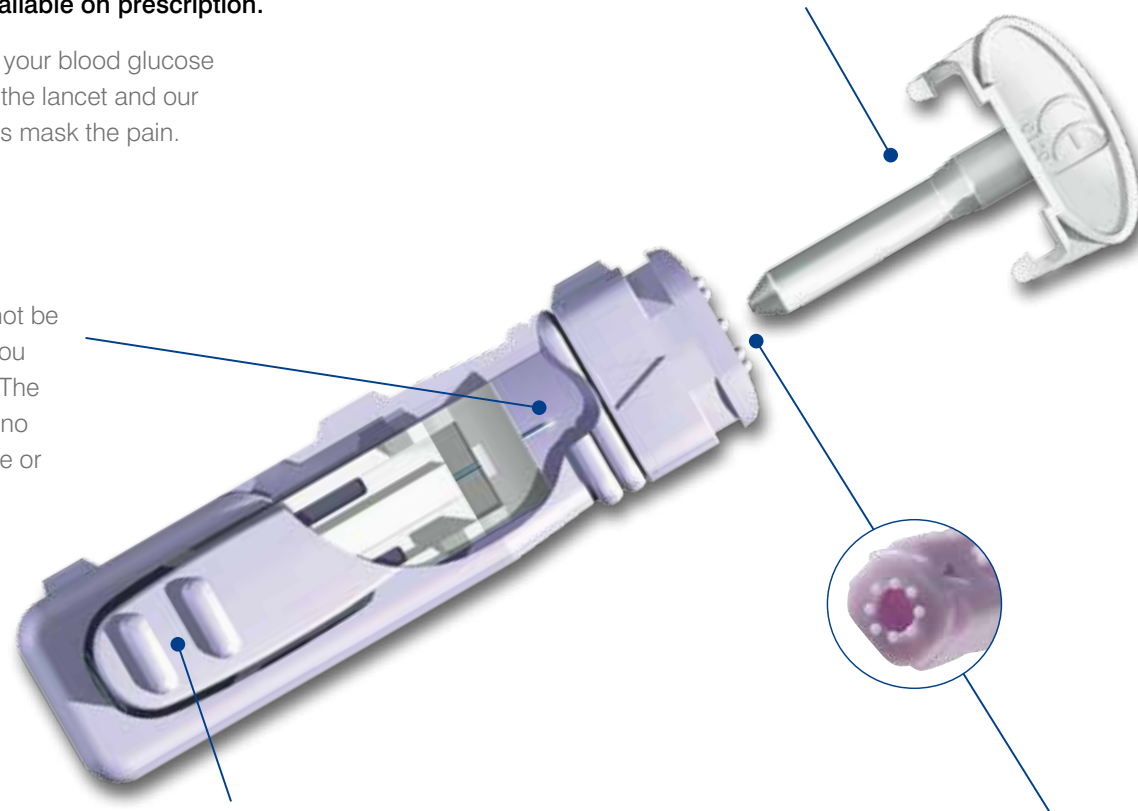

Simple.

One-click operation.


Comfort Zone Technology® takes your mind of the pain.

Pressing Unistik® 3's 'Comfort Zone' of 8 raised dots against the skin sends a message of comfort to the brain. Then, when you click Unistik® 3's release button, the lancet penetrates your skin, but the dominant Comfort Zone sensation means the pain 'message' to your brain is greatly reduced. Distracting the brain in this way makes sampling virtually pain-free.




How to use Unistik® 3


Using Unistik® 3 couldn't be simpler

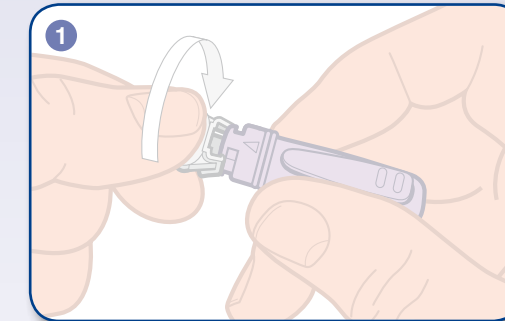

TWIST off and discard lancet cap.

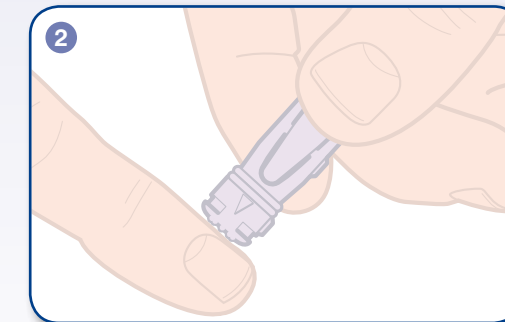

Press Unistik® 3 against finger and **PRESS** the release button.

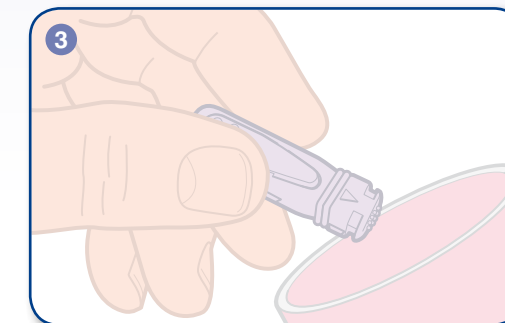

DISPOSE of the used lancet in a suitable sharps container.

For further information please call us on 0800 731 6959.